



All-In-One Device



Core Body
Thermometer



Reflective Pulse
Oximeter



Digital
Stethoscope



Single Lead
ECG



Ultra HD
Camera



Third Party
Sensors*

*Third party sensor data such as Blood Pressure, Glucometer, Spirometer, Weight, 12-Lead ECG and Ultra Sound.



www.medwandhealth.com/contact

Clinical Care, Anywhere

Whether in an On-Site exam, or a Remote Consultation, MedWand enhances clinical care by providing live monitoring, remote clinical-grade vitals capture, and support for external AI-based vitals data analysis.

TYPICAL USE CASES



Nursing



Rural Clinic



Corporate



Transport



Education



Clinical Trials

Integrated Workflow



01.

Launches
from Your
HIS/EHR




02.


Data Captured
by VirtualCare
Software




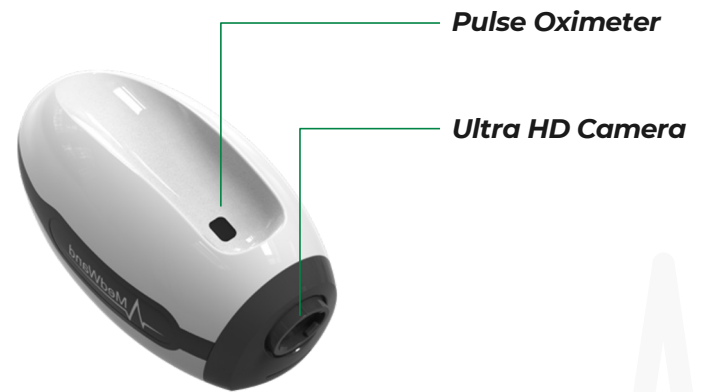
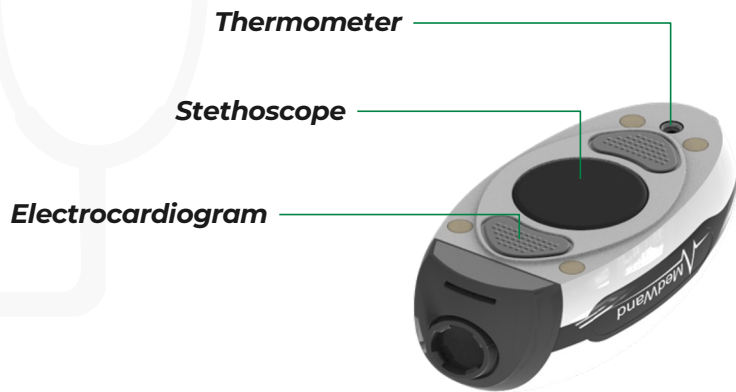
03.

Review and
Export Data
to HIS/EHR

 NORTH LAS VEGAS, NV

 info@medwand.com

 +1-844-MED-WAND



8+
Medical Tools



Ear / Nose



Throat



Skin

In-Field Application can also Ingest 3rd Party Device Data such as:

Blood Pressure


Glucometer

Spirometer

Weight

12-Lead ECG


MEDWAND EXAM DATA WORKFLOW

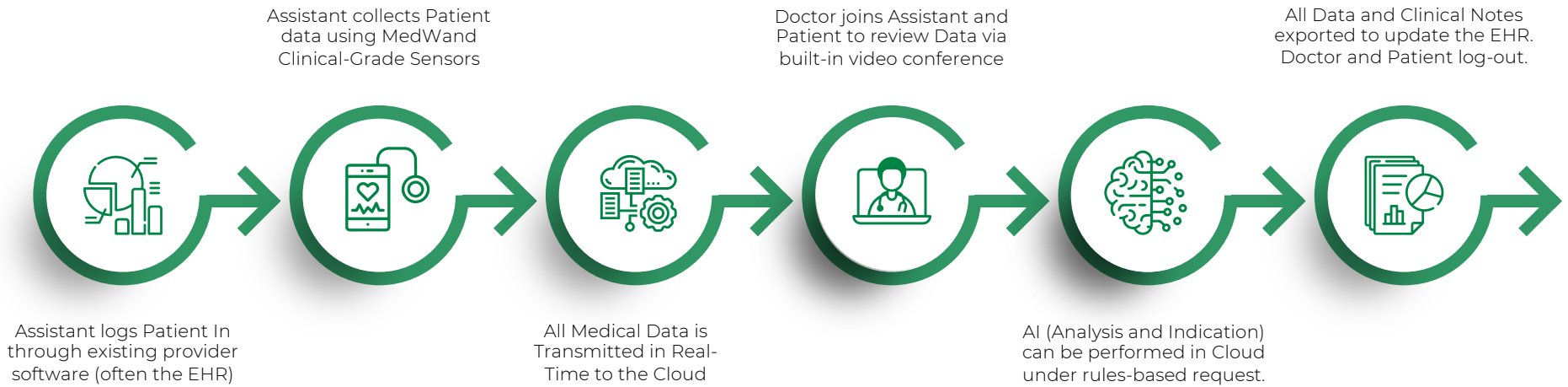


MEDWAND EXAM INTERFACE EXAMPLES



Designed to operate in On-Site (local), Virtual Consultation (remote) or Hybrid (start local then remote) exam modes.

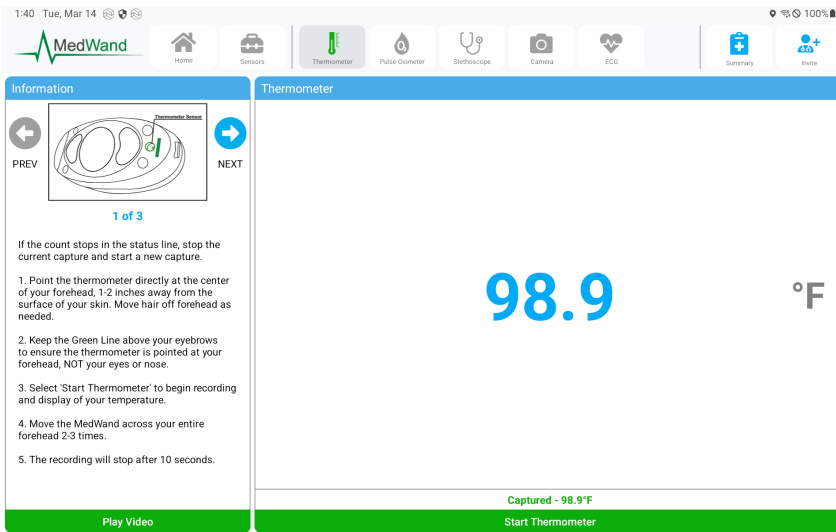

Patient/Nurse view during On-Site local exam mode

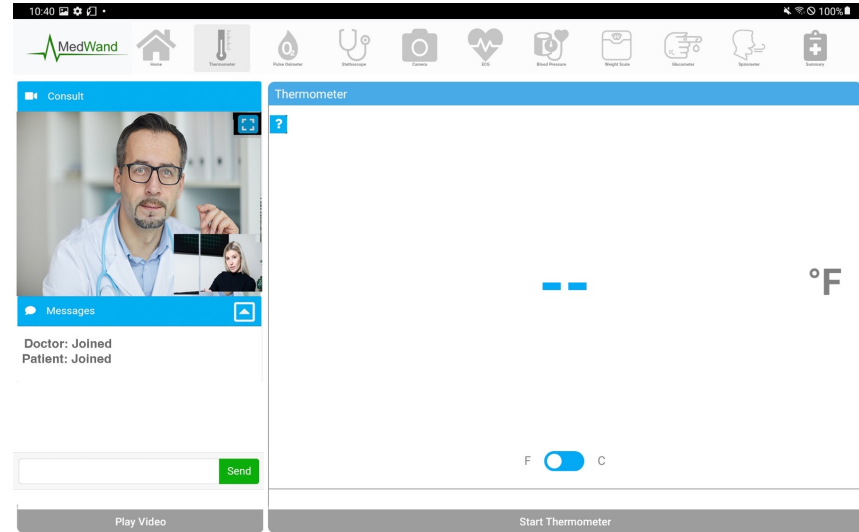

Patient/Nurse view during Virtual Consult remote exam mode

In Hybrid mode, Patient/Nurse uses Invite icon to convert an existing On-Site exam into a Virtual Consultation with supervising clinician.

MEDWAND CLINIC CONFIGURATIONS



HOME CLINIC

Designed for concierge, work, and international travel applications.



MOBILE CLINIC

Designed for home visits, assisted living, and skilled nursing facilities.



RUGGEDIZED CLINIC

Designed for rural health care, air or sea transportation, and industry.



DESKTOP CLINIC


Designed for corporate wellness, pharmacies, and office-based clinics.

MEDWAND KIOSK SERVICES



Supports Exam Appointments, Kiosk-style Exam Access, Exam Result Review, Sharing Exam as HTML file.

VirtualCare Kiosk : Appointments				
Code	Type	Patient Name	Clinician Name	Date Created ^
138755	Remote Exam	Alger, Bianca (96059)		
270377	OnSite Exam	Alger, Bianca (96059)		
406419	Patient Exam	Alger, Bianca (96059)		
944690	OnSite Exam	Doe, Jane		
638883	Remote Exam	Alger, Bianca (96059)		

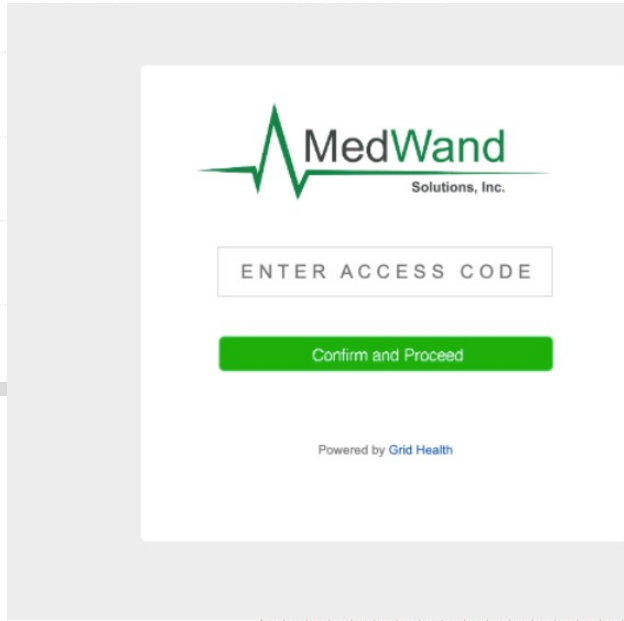

MedWand Solutions, Inc.

ENTER ACCESS CODE

Confirm and Proceed

Powered by Grid Health



7/15/2023 12:22:16 AM


Paul Simmonds 100 5634 223

Practitioner: Jane Porter, NP

Sensor Readings

PainLevel 5 Numerical	PulseRate 74	Respiratory 15.0	SpO2 95	Thermometer 96.9 F
-----------------------------	-----------------	---------------------	------------	-----------------------

-14 17:19:58 UTC-07:00

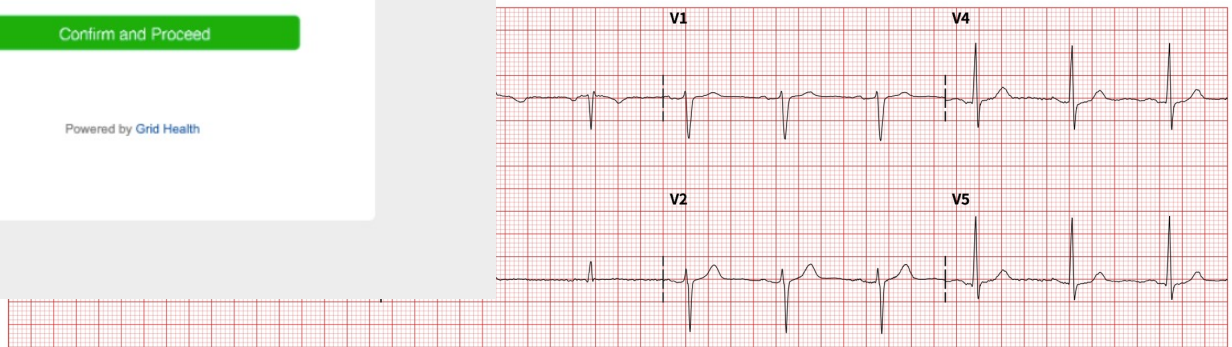






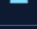



MEDWAND COMMAND CENTER










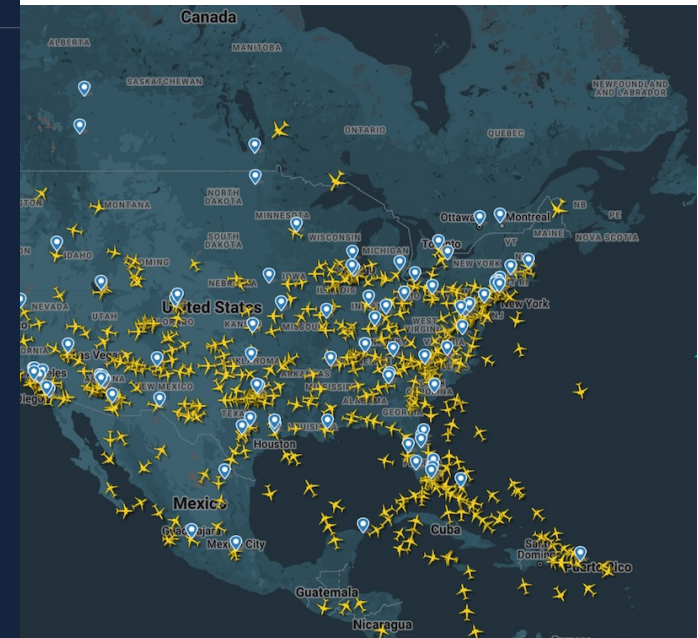
Supports Real-Time Exam Monitoring, provides Response Center for Exam Escalation Requests

MedWand Command Center

Search...

-  Martha Edw, Female, 61
TC-FDC > TK-73 1
-  Thomas Yoo, Male, 117
TC-FDC > TK-48 1
-  Shelia Glo, Female, 45
TC-FDC > TK-77 1
-  Daniel Wal, Male, 56
TC-FDC > TK-57 1
-  Stephanie, Female, 4
TC-FDC > TK-63 1
-  Howard Ben, Male, 32
TC-FDC > TK-15 1
-  William Pe, Male, 0
TC-FDC > TK-4 1
-  Rachael Pa, Female, 49
TC-FDC > TK-11 1
-  Chrystal R, Female, 16 1

 Lyn Glover, Female, 64 TC-FDC > TK-1 2 <table border="1" style="width: 100%; border-collapse: collapse; font-size: 8px;"> <tr><td>HR</td><td>BP</td><td>SpO2</td></tr> <tr><td>83</td><td>114/87</td><td>89</td></tr> <tr><td>T</td><td>RR</td><td>Glc</td></tr> <tr><td>37 C</td><td>14</td><td>194</td></tr> <tr><td>PEF</td><td>FEV1</td><td>FVC</td></tr> <tr><td>16</td><td>94</td><td>16</td></tr> <tr><td>Stethoscope</td><td>Othoscope</td><td>ECG</td></tr> </table>	HR	BP	SpO2	83	114/87	89	T	RR	Glc	37 C	14	194	PEF	FEV1	FVC	16	94	16	Stethoscope	Othoscope	ECG	 Charles Kl, Male, 30 TC-FDC > TK-2 2 <table border="1" style="width: 100%; border-collapse: collapse; font-size: 8px;"> <tr><td>HR</td><td>BP</td><td>SpO2</td></tr> <tr><td>83</td><td>114/87</td><td>89</td></tr> <tr><td>T</td><td>RR</td><td>Glc</td></tr> <tr><td>37 C</td><td>14</td><td>194</td></tr> <tr><td>PEF</td><td>FEV1</td><td>FVC</td></tr> <tr><td>16</td><td>94</td><td>16</td></tr> <tr><td>Stethoscope</td><td>Othoscope</td><td>ECG</td></tr> </table>	HR	BP	SpO2	83	114/87	89	T	RR	Glc	37 C	14	194	PEF	FEV1	FVC	16	94	16	Stethoscope	Othoscope	ECG	 Bertha Tom, Female, 13 TC-FDC > TK-3 2 <table border="1" style="width: 100%; border-collapse: collapse; font-size: 8px;"> <tr><td>HR</td><td>BP</td><td>SpO2</td></tr> <tr><td>123</td><td>113/37</td><td>92</td></tr> <tr><td>T</td><td>RR</td><td>Glc</td></tr> <tr><td>35 C</td><td>10</td><td>216</td></tr> <tr><td>PEF</td><td>FEV1</td><td>FVC</td></tr> <tr><td>91</td><td>84</td><td>47</td></tr> <tr><td>Stethoscope</td><td>Othoscope</td><td>ECG</td></tr> </table>	HR	BP	SpO2	123	113/37	92	T	RR	Glc	35 C	10	216	PEF	FEV1	FVC	91	84	47	Stethoscope	Othoscope	ECG	 William Pe, Male, 0 TC-FDC > TK-4 1 <table border="1" style="width: 100%; border-collapse: collapse; font-size: 8px;"> <tr><td>HR</td><td>BP</td><td>SpO2</td></tr> <tr><td>116</td><td>115/80</td><td>99</td></tr> <tr><td>T</td><td>RR</td><td>Glc</td></tr> <tr><td>37 C</td><td>13</td><td>301</td></tr> <tr><td>PEF</td><td>FEV1</td><td>FVC</td></tr> <tr><td>63</td><td>33</td><td>6</td></tr> <tr><td>Stethoscope</td><td>Othoscope</td><td>ECG</td></tr> </table>	HR	BP	SpO2	116	115/80	99	T	RR	Glc	37 C	13	301	PEF	FEV1	FVC	63	33	6	Stethoscope	Othoscope	ECG
HR	BP	SpO2																																																																																					
83	114/87	89																																																																																					
T	RR	Glc																																																																																					
37 C	14	194																																																																																					
PEF	FEV1	FVC																																																																																					
16	94	16																																																																																					
Stethoscope	Othoscope	ECG																																																																																					
HR	BP	SpO2																																																																																					
83	114/87	89																																																																																					
T	RR	Glc																																																																																					
37 C	14	194																																																																																					
PEF	FEV1	FVC																																																																																					
16	94	16																																																																																					
Stethoscope	Othoscope	ECG																																																																																					
HR	BP	SpO2																																																																																					
123	113/37	92																																																																																					
T	RR	Glc																																																																																					
35 C	10	216																																																																																					
PEF	FEV1	FVC																																																																																					
91	84	47																																																																																					
Stethoscope	Othoscope	ECG																																																																																					
HR	BP	SpO2																																																																																					
116	115/80	99																																																																																					
T	RR	Glc																																																																																					
37 C	13	301																																																																																					
PEF	FEV1	FVC																																																																																					
63	33	6																																																																																					
Stethoscope	Othoscope	ECG																																																																																					
 Katherine, Female, 99 TC-FDC > TK-5 2 <table border="1" style="width: 100%; border-collapse: collapse; font-size: 8px;"> <tr><td>HR</td><td>BP</td><td>SpO2</td></tr> <tr><td>68</td><td>113/55</td><td>95</td></tr> <tr><td>T</td><td>RR</td><td>Glc</td></tr> <tr><td>37 C</td><td>10</td><td>162</td></tr> <tr><td>PEF</td><td>FEV1</td><td>FVC</td></tr> <tr><td>49</td><td>66</td><td>95</td></tr> <tr><td>Stethoscope</td><td>Othoscope</td><td>ECG</td></tr> </table>	HR	BP	SpO2	68	113/55	95	T	RR	Glc	37 C	10	162	PEF	FEV1	FVC	49	66	95	Stethoscope	Othoscope	ECG	 Aaron Sals, Male, 55 TC-FDC > TK-6 2 <table border="1" style="width: 100%; border-collapse: collapse; font-size: 8px;"> <tr><td>HR</td><td>BP</td><td>SpO2</td></tr> <tr><td>102</td><td>110/64</td><td>97</td></tr> <tr><td>T</td><td>RR</td><td>Glc</td></tr> <tr><td>37 C</td><td>7</td><td>73</td></tr> <tr><td>PEF</td><td>FEV1</td><td>FVC</td></tr> <tr><td>42</td><td>80</td><td>47</td></tr> <tr><td>Stethoscope</td><td>Othoscope</td><td>ECG</td></tr> </table>	HR	BP	SpO2	102	110/64	97	T	RR	Glc	37 C	7	73	PEF	FEV1	FVC	42	80	47	Stethoscope	Othoscope	ECG	<p>TC-FDC > TK-7 Hypertensive Crisis!</p>	<p>TC-FDC > TK-8 High Glucose Alert!</p>																																										
HR	BP	SpO2																																																																																					
68	113/55	95																																																																																					
T	RR	Glc																																																																																					
37 C	10	162																																																																																					
PEF	FEV1	FVC																																																																																					
49	66	95																																																																																					
Stethoscope	Othoscope	ECG																																																																																					
HR	BP	SpO2																																																																																					
102	110/64	97																																																																																					
T	RR	Glc																																																																																					
37 C	7	73																																																																																					
PEF	FEV1	FVC																																																																																					
42	80	47																																																																																					
Stethoscope	Othoscope	ECG																																																																																					
<p>TC-FDC > TK-9 High Heart Rate!</p>	<p>TC-FDC > TK-10 Low Oxygen Saturation!</p>	<p>TC-FDC > TK-11 Critical - Hypothermia!</p>	 Susanne Sm, Female, 119 TC-FDC > TK-12 2 <table border="1" style="width: 100%; border-collapse: collapse; font-size: 8px;"> <tr><td>HR</td><td>BP</td><td>SpO2</td></tr> <tr><td>102</td><td>114/39</td><td>85</td></tr> <tr><td>T</td><td>RR</td><td>Glc</td></tr> <tr><td>38 C</td><td>8</td><td>341</td></tr> <tr><td>PEF</td><td>FEV1</td><td>FVC</td></tr> <tr><td>25</td><td>88</td><td>91</td></tr> <tr><td>Stethoscope</td><td>Othoscope</td><td>ECG</td></tr> </table>	HR	BP	SpO2	102	114/39	85	T	RR	Glc	38 C	8	341	PEF	FEV1	FVC	25	88	91	Stethoscope	Othoscope	ECG																																																															
HR	BP	SpO2																																																																																					
102	114/39	85																																																																																					
T	RR	Glc																																																																																					
38 C	8	341																																																																																					
PEF	FEV1	FVC																																																																																					
25	88	91																																																																																					
Stethoscope	Othoscope	ECG																																																																																					



Map showing patient locations across North America (Canada, United States, Mexico, Central America). Yellow airplane icons indicate patient locations.

MEDWAND CONTACT INFORMATION



www.medwand.com

1-844-MEDWAND

info@medwand.com

MedWand™ Specifications

General:

- Intended for use by adults (18yrs or older).
- Non-cytotoxic for all patients and user contact per ISO 10993-5.
- Non-irritating and non-sensitizing for all patient and user contact per ISO 10993-10.
- Power input USB type C at 5 volts DC, 500 mA max current, 2.5 watts max power.
- Ingress protection rating of IP22 per IEC 60529-1 and IEC 60601-1-11.
- All materials compliant to RoHs2 directive 2011/65/EU.
- All materials compliant with WEEE directive 2012/19/EU to support e-waste safety.
- Weight 5.75oz, Height 2.40in, Width 2.5in, Length 5.00in (163g, 61 x 64 x 127mm)

Environmental:

- Operating Temperature: 61°F to 104°F (16°C to 40°C).
- Operating Relative Humidity: 20% to 95% non-condensing.
- Operating Altitude: Sea level to at least 3,000 meters.
- Storage Temperature: - 4°F to 122°F (-20°C to 50°C).
- Storage Relative Humidity: 5% to 95% non-condensing.

Pulse Oximeter:

- Reflective-mode Pulse Oximeter reporting SpO2 and Pulse Rate values.
- Robust resilience to fingertip movement providing high Signal to Noise ratio.
- Accuracy meets requirements of ISO 80601-2-61:2017 or later version between 75% and 100% blood oxygen saturation.



Electrocardiogram:

- Single lead dual pad touch only system.
- Equivalent to 12 lead ECG Lead I (RA/LA).
- Less than 2.2% total gain error.
- Accuracy meets requirements of IEC 60601-2-47:2011 or later version.
- Recorded in MedWand app as a BMP file.



IR Thermometer:

- Non-contact infrared sensor.
- Lab tested accurate to +/- 0.3°C.
- Projects accurate core temperature. between 94°F and 107°F (34°C and 44°C).
- Operates at ambient temperatures between 61°F and 104°F (16°C and 40°C).
- Accuracy meets requirements of ISO 80601-2-56:2017 or later version.



Stethoscope:

- Digital audio output at 44.1kHz sample rate.
- Frequency response 5Hz - 15kHz.
- User selectable osculation audio filters;
 - Heart (low pass)
 - Lungs (high pass)
 - Abdomen (wide band pass)
- Lab verified performance equal or superior to most available digital stethoscopes.
- Recorded in MedWand app as a WAV file.



Ultra HD Camera:

- 5.7mm CMOS sensor.
- 3200x 2400 max resolution
- Auto and manual focus modes.
- Field of view 55.4°.
- Total distortion <1%.
- User adjustable LED illumination.
- Includes adaptors and tongue depressor slot for Otoscope, Dermatoscope and throat examination still image capture functions.
- Accepts Welch Allyn 52134, 52133, 52135 Specula.
- Accepts standard wood tongue depressor.
- Recorded in MedWand app as JPEG file.





Specifications:

Length - 5 inches


Width - 2.4 inches

Height - 2.5 inches

Weight - 5.750 ounces

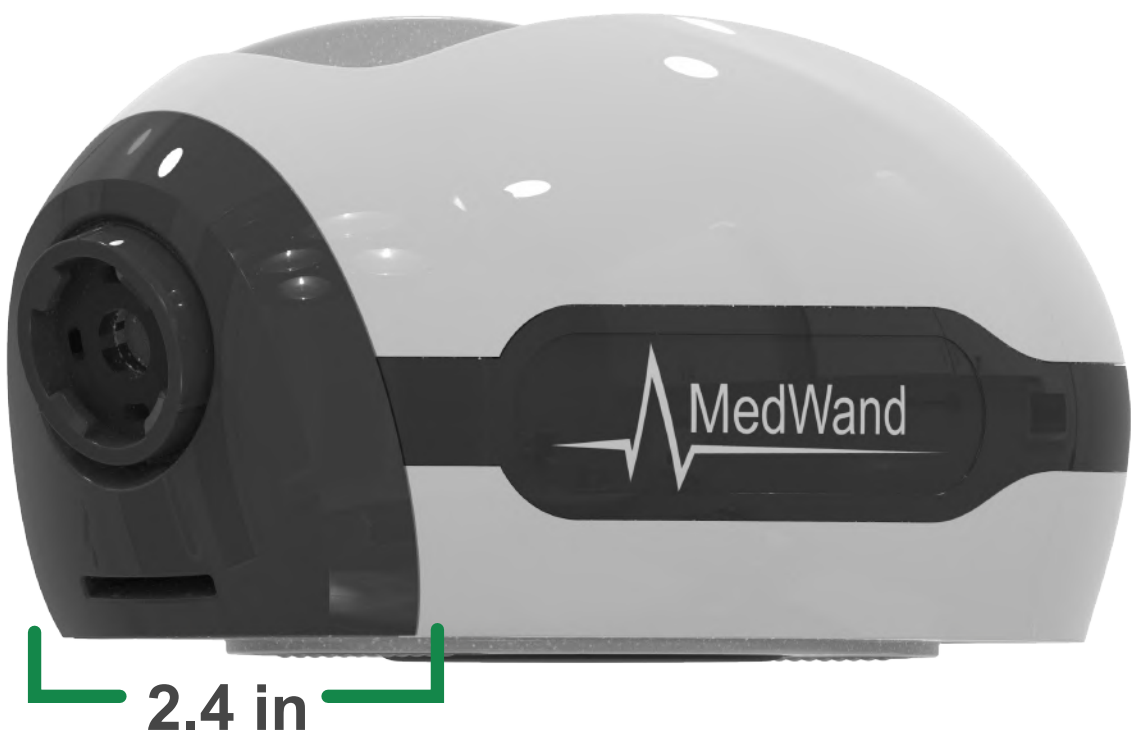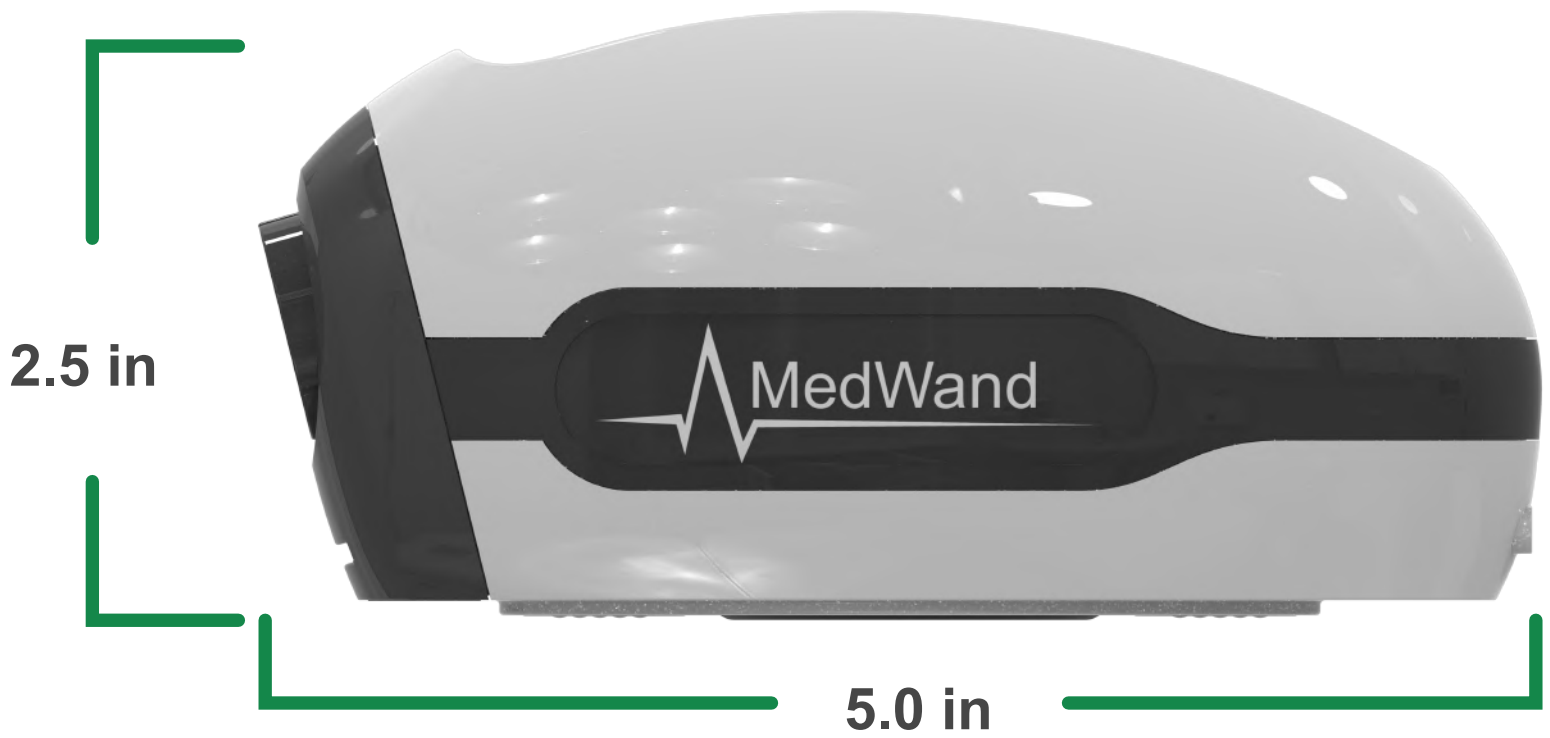

DESKTOP CLINIC

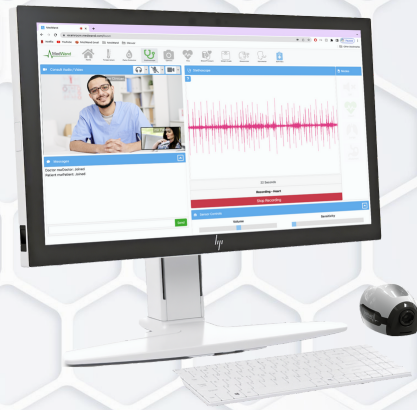

Includes:

VirtualCare software
10" Android tablet
Desktop UV sanitizer

Optional Blood Pressure
monitor, Spirometer and
Glucometer

Healthcare for the 21st Century





Command Center Option

Enable one care center clinician to monitor hundreds of exam results.

Flexible Exam Escalation

Convert On-Site exam into a full Remote Consultation as needed.

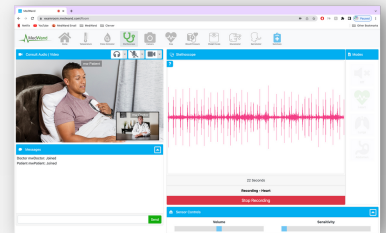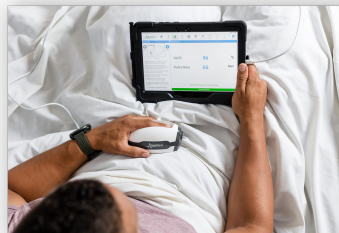

📍 NORTH LAS VEGAS, NV

✉ info@medwand.com

☎ +1-844-MED-WAND